







LED color light skin beauty

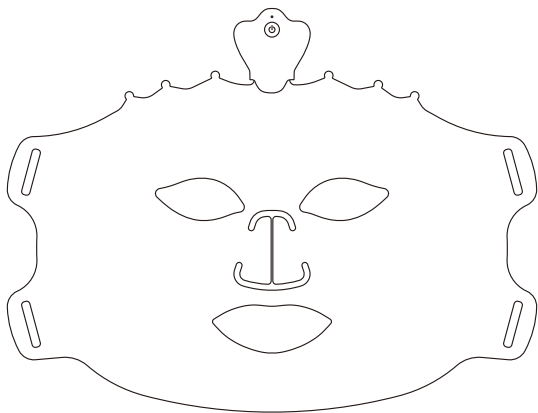

Initial setup

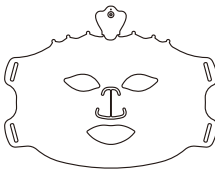
Before opening the package for use, be sure to charge the handle of the mask controller fully. The indicator light of the controller will start to flash, indicating that it is being charged. After the battery is fully charged, the battery indicator stops blinking. The charging time is 4-5 hours. The adjustable velcro allows you to wear the mask comfortably. To install the mask, pass the strap through the entrance hole on the other side of the mask, wrap the strap around your head and secure it firmly.

Product list

1 	2 	3 
4 	5 	1. Beauty device x1 2. Eye mask x2 3. Charging cable x1 4. Fasten strap x2 5. Instructions x1

PRODUCT EFFICACY

LED color light beauty facial mask instrument is based on the principle of light irradiation, safe and simple, and can be reused to maintain the facial skin. It uses a technique of quadruple composite phototherapy, combining red, yellow, blue, and far-infrared light into three main modes. You can freely choose the modes, including Acne removing mode, soothing mode, and Whitening and brightening mode. Four fold light wave coordinated repair, with overlapping effects and double energy.

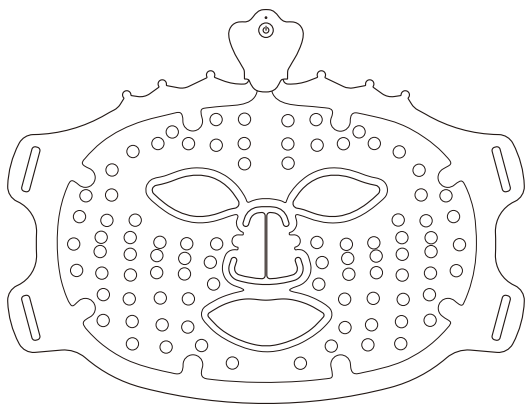

LIGHTING INTRODUCTION

Red light (wavelength 630nm):

Known as "bioactive light". It can improve cell activity, accelerate metabolism and promote collagen formation, with whitening and light spots, skin rejuvenation and wrinkle removal, increase skin elasticity. improve skin condition, repair damaged skin, reduce pores and tighten skin, anti-aging and antioxidant effects.

Blue light (wavelength 470nm):

It can efficiently destroy Propionibacterium and other bacteria with no damage to the skin tissue, and reduce the inflammation of acne to heal in a short time. It is effective in removing acne, sterilization and anti-inflammation, repairing skin without leaving acne marks and pigmentation.

Green light (wavelength 520nm):

Neutralizing and balancing, it can reduce skin oil secretion, balance the ratio of water to oil, effectively relieve mental tension, lymphatic drainage and edema.

Yellow light (590nm wavelength after special treatment):

Improve the exchange of cellular oxygen, give skin cells to replenish energy, decompose pigment, promote lymphatic detoxification. improve skin roughness, wrinkles, treat skin redness and fever, moss etc, enhance immunity.

Purple light:

It combines the efficacy of two kinds of light therapy, especially in the treatment of acne and acne scars, and has a particularly good effect and repair effect.

Green and blue light:

It can gradually enhance cellular energy and has a very good effect on metabolism.

Laser light (white light):

Penetrates deep into the skin, accelerates the metabolism of active tissue, decomposes pigmentation, improves fine lines and sagging skin.

Instructions

Step 1:

Clean and dry your face.



Step 2:

Attach the adjustment strap, wear the mask, adjust the strap strap to a comfortable position.



Step 3:

Long press the controller button for 1.5 seconds to turn on the photon.
Tap to switch, freely choose your skin care mode.

Switching photon switches
Long press to power on, beep once
Long press to shut down, beep twice



Step 4:

After 10 minutes, the device automatically shuts down, remove the mask and follow up with skin care.
(Please moisturize your skin immediately after use)



Before use, please first irradiate the face for 1 minute, and continue to use without discomfort. Please start from the low level to establish skin tolerance and choose the energy gear according to the actual situation of the skin.

Do not increase the illumination time and frequency blindly, please adjust the illumination time according to the actual condition of the skin.

It is recommended that the frequency of use of this product is 3-5 times a week, and the single use time is 10 minutes , 4 consecutive weeks a cycle.

When using for the first time, make sure the device is fully charged.

If the indicator blinks, the charging is ongoing. If the indicator blinks, the charging is complete.



Do not use the apparatus during charging.

CAUTIONS:

1. Please use the product strictly in accordance with the instructions, do not use the apparatus beyond the instructions, otherwise it may lead to instrument failure or accidents.
2. Do not use this product to improve any condition other than the indication.
3. Pregnant women, breastfeeding and children prohibited from using this product.
4. If the user suffers from physical, sensory and mental illness, or is undergoing dermatological improvement including the onset period (such as lupus, photosensitive eczema, or albinism), this product is not recommended to use, and it should only be used under supervision or under the guidance of a responsible person (such as a dermatologist).
5. Do not use this product if you suffer from any photosensitivity disorder (photosensitivity), otherwise, you may have a severe skin reaction.
6. Do not use on the same day as skin care products containing photosensitive ingredients.
7. Pay attention to eye protection during use, put on the eye mask correctly, it is recommended to close the eyes during use.
8. If there is any discomfort during the use of this product, please stop using it immediately.
9. This product is not waterproof, do not immerse the equipment in water or wash with water.
If it has been in contact with water, stop use it.
10. Do not use for infants and children, and store out of reach of children.
11. If any parts of the instrument are damaged, please do not use this product.
12. This product is for daily use only. Please do not disassemble, transform or repair it without authorization.
13. It is recommended that the use rate of this product is 3-4 times a week, and the single use time is controlled at about 10 minutes.
14. Please timely moisturize the skin after use.

PRODUCT PARAMETERS :

Name	LED color light skin beauty Contour face
Number of lamp beads	80x3=240 beams
Time timing	10min
Adapter Input	220V
Outer box size	338x258x45mm
Net weight	240g
Battery capacity	500mAh
Adapter output	5V1A



