

Caution

1. To move this instrument, it is necessary to lift it up and put it down gently; in case of the failure, please do not knock hard on it so as to prevent malfunction caused by the impact.
2. To move this instrument, it is necessary to keep it horizontal so as to avoid vibration and impact.
3. Do not squeeze the instrument or place heavy stuff onto the instrument so as to prevent shell deformation, damage or malfunction.
4. It is necessary to select the corresponding working power supply in accordance with the marks on the instrument body.
5. It is not allowed to use this instrument in the bathroom.
6. To prevent the lightning accident, when it thunders, please show down the power switch and pull out the power plug.
7. It is not allowed to inject any fluid into the instrument body so as to prevent the instrument from being damaged due to the short circuit.
8. If the instrument will not be used for a long time, it is necessary to pull out the power plug from the socket.

Technical parameters

1. Operating voltage: AC 220V--AC 240V
2. Maximum power of the complete machine: 20W
3. Packing size: 27*14*27 (L*W*H) cm

Version: 20140822

Microcomputer Micro Dermabrasion Beauty Instrument User Manual




Dear customers,

Thank you for choosing our product! Please go through this manual before use in order to ensure the maximum performance of the instrument and the safety in use. Please keep properly for further requires.

Function

Get rid of aging skin, improve skin texture, whitening, improve scars and skin irregularities; accelerate blood circulation and improve large pores

Operation Panel Description

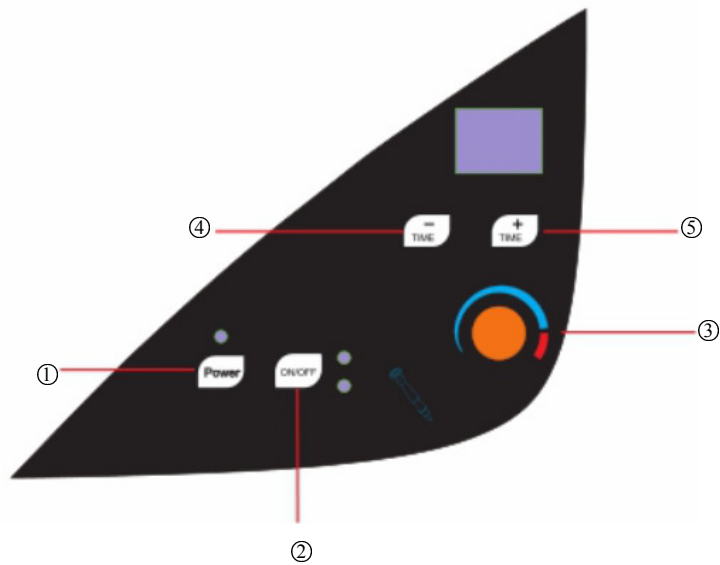

- ① Power ② ON/OFF ③ Pressure Regulator
④ Time Decrease ⑤ Time Increase

Operation

1. Choose the suitable tip and install on the top of the handle, connect the handle into the interface behind the machine with hose.
2. Plug in, press ①, indicator turns green, time shows 30 (minutes), and the indicator of ② turns red, which means the machine stands by.
3. Press ②, the indicator turns green and the second point of the time display flashes, which means the machine is starting working.
4. Choose the suitable strength from low to high by adjusting the pressure regulator in clockwise direction.

Notes

1. If beautician does the treatment with sterilized gloves, hands should first be cleaned and disinfected with alcohol
 2. The accessories should be sterilized and keep clean and dry
- The machine is equipped with dermabrasion tips of various sizes for different part use. Installation as below photo:

