
EN

Instructions Manual

Podiatry Bed

2232A.1.A26.N

Instruction Manual

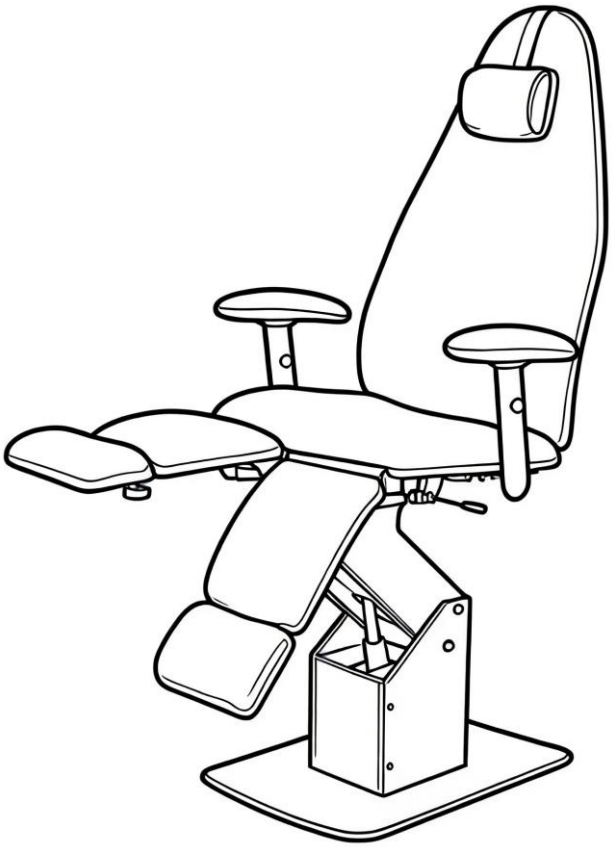VERSION 1.0

Thank you very much for choosing our product. This manual is compiled to help operators use the product properly. We hope our product and service bring you a satisfying experience.

This user manual contains general information about this product. It serves as a starting point for new users, giving a general introduction to the system, its working principle, installation procedures and various operation methods.

Important Warnings

1. The power supply must be disconnected when replacing electrical components, maintaining or cleaning the unit. The power plug shall be connected to a properly installed and earthed socket to prevent fire hazard and electric shock risks.
2. Do not stand on this equipment.
3. Do not allow children to play with the equipment. Keep the unit out of children's reach. Adjust the backrest to its highest position before leaving.
4. Do not place the equipment in wet locations such as bathrooms to avoid electric shock or equipment failure. Keep controllers away from water contact.
5. Do not place or operate the equipment in areas where painting/coating work is in progress or where oxygen is being supplied.
6. Persons with physical impairments, tactile disorders, mental disorders, lack of relevant experience or common sense, as well as children, shall NOT use this equipment unless supervised or instructed by personnel who take responsibility for their personal safety.
7. If the power plug is damaged, contact your distributor, service centre or qualified technical personnel for repair to prevent hazards.
8. Ensure no obstacles are present within the full movement stroke during any operation. The equipment shall only be operated by trained professionals.

Precautions

1. The equipment shall be managed by experienced personnel who are familiar with the product's structure, performance and operating procedures. Regular maintenance shall be performed.
2. Read this manual thoroughly before installation, operation or adjustment, and keep it in a safe place. Do not apply load to only one side or single point of this product.
3. Follow safety operating rules. Do not overload the product. Make sure all safety protection components are complete and functional, and eliminate potential safety hazards in a timely manner.
4. Tighten bolts and nuts on a regular basis to prevent malfunctions caused by loose, slipping, vibrating or falling off parts. Do not arbitrarily disassemble any bed components during maintenance work.
5. Never operate the equipment when a fault occurs. Stop operation immediately if you notice abnormal phenomena or unusual noises. Do not resume use until the fault has been located and resolved.
6. Operators must be familiar with all precautions stated in this manual to guarantee personal and equipment safety.
7. This electric podiatry bed is fitted with a 3 core power cord. It must be plugged into a properly installed 3 pin earthed socket. This prevents electric shock and ensures proper bed operation.

8. Operate the equipment strictly in accordance with warning labels affixed to the product.
9. When controlling the full bed lifting movement, keep holding the button continuously until the bed reaches your desired position.
10. Motor duty cycle limit: Maximum continuous ON time: 2 minutes; minimum OFF time: 18 minutes.
11. All maintenance work shall be carried out by trained professionals only.
12. After each treatment, wipe dust and stains off leather and wooden components using a neutral cleaning agent.
13. When adjusting the backrest, seat cushion, footrest or overall bed height via hand controller or optional foot controller, keep your hands, feet and head clear of all moving parts.

Scope of Application

This product is suitable for beauty treatment, massage and pedicure / manicure services.

Installation

1. Position the bed according to the layout of your salon or home, lighting conditions and practical working requirements. Keep it in a clean, dry and cool environment.
2. Unpack the carton. Inspect the bed for physical damage and verify all accessories and spare parts are complete. Contact your distributor or manufacturer if any issues are found.
3. Place the bed on firm, level ground to prevent accidents.
4. Connect the bed's power plug correctly to the mains power supply.

Product Structure

This equipment consists of a control box, electric motor, metal hardware, plastic parts and upholstery leather. The control box operates motor functions and communicates with the hand-controller. This integrated unit supports manual backrest-angle adjustment, manual footrest-angle adjustment, and motor-driven full-bed height adjustment.

Controllers for Model 2232A.1.A26.N

Standard Configuration: Hand-Controller

No	Function	Operation
1	Bed move upward	Press and hold the button; the bed rises smoothly
2	Bed move downward	Press and hold the button; the bed lowers smoothly

Maintenance

1. After every treatment session, clean dust and stains on leather upholstery, metal parts and plastic covers only with neutral cleaning solutions. Do NOT use alcohol-based, acidic or alkaline disinfectants, as these will shorten the service life of leather and plastic components.
2. Lubricate all moving mechanical parts once per year.

Circuit System Diagram

Power Supply → Control Box Hand Controller (Standard) / Foot Controller (Optional) connect to Control Box Control Box drives Full bed Lifting Motor


Troubleshooting

Fault Phenomenon	Possible Causes	Solutions
Bed cannot lift or lower	1. Poor power plug connection4. Defective rocker switch5. Faulty control box6. Broken motor tail wiring / connector7. Damaged motor	1. Re connect and firmly insert the power cord / plug.2. Reposition movable footrest components away from button interference; open fixed footrest housing, inspect and re secure power wiring terminals.3. If the problem persists after above steps, contact your service centre or supplier.
Abnormal operating noise	1. Insufficient lubrication on rotating components2. Insufficient lubrication for POM plastic assemblies3. Motor quality defect	1 2: Apply suitable lubricant to relevant parts. For motor related defects, contact after sales support.


If faults remain unresolved after performing the above checks, please contact your service centre or supplier for assistance.

Diagram Reference List


1. 2232A (1 motor) Main unit structural drawing ① Backrest, manual angle adjustment ② Full bed motor driven lifting ③ Footrest, manual angle adjustment


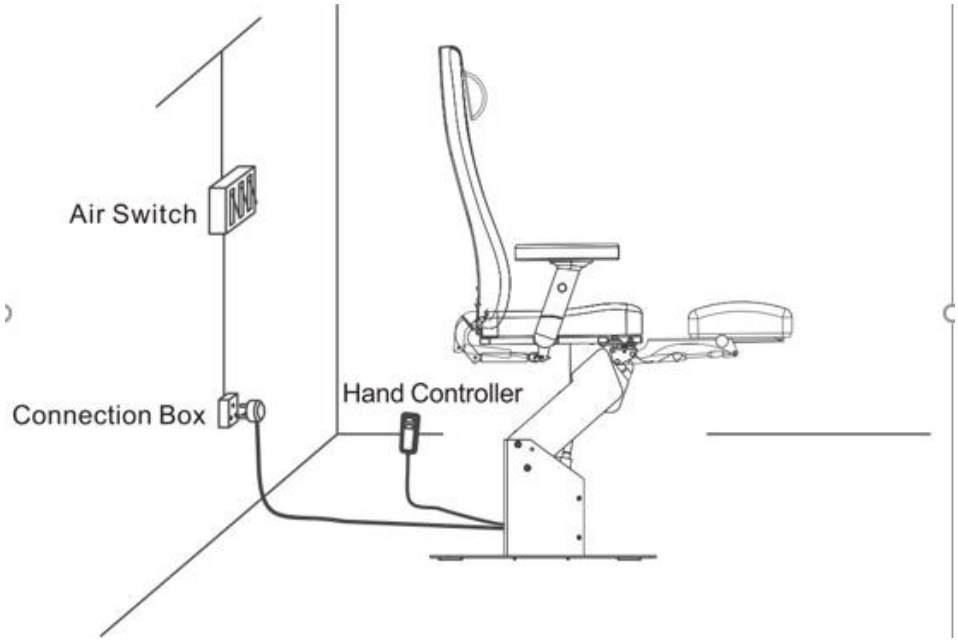
2. Whole unit installation & wiring schematic Includes: Air Switch, Connection Box, Hand Controller, Podiatry Bed main body

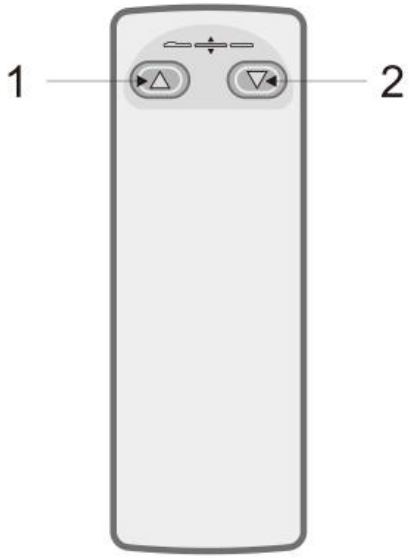
3. Standard Hand controller drawing 1 = UP button, 2 = DOWN button; buttons must be held continuously during adjustment

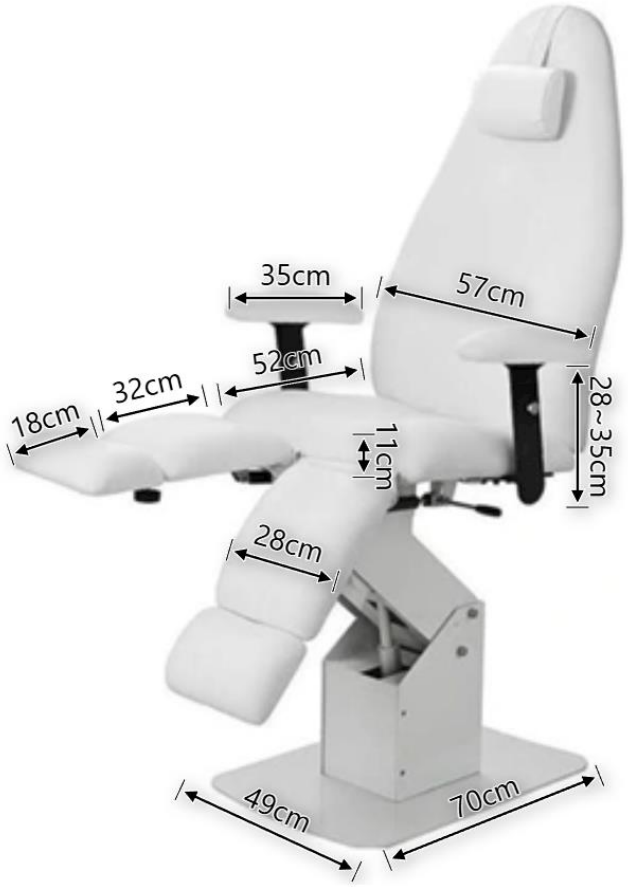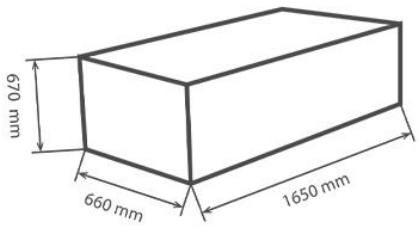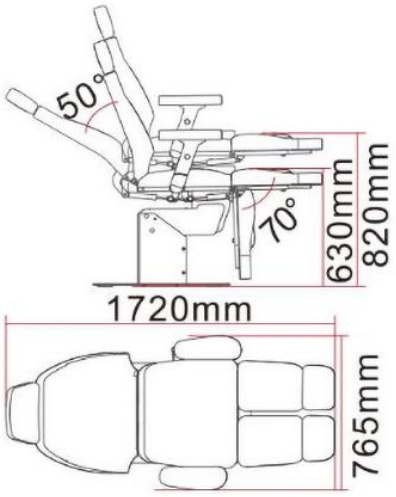4. Packing size & dimension drawing for 2232 Showing carton parameters, net/gross weight, bed outline dimensions, backrest & footrest rotation angles

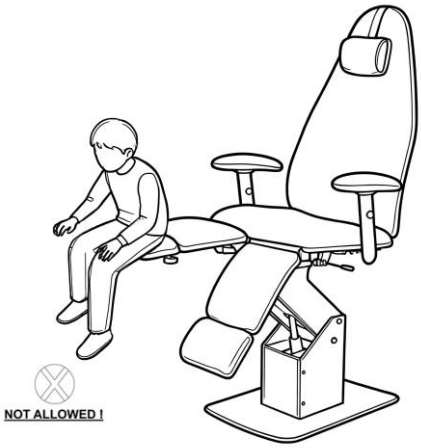5. Safety prohibition drawing ⚠ NOT ALLOWED: Sitting or applying load only on the far side edge of the bed.

Technical Specifications

Item	Specification
Power Supply	AC 100 120V / AC 220 240V / AC 100 240V, 50/60Hz
Max. Power Consumption	300 W
Ambient Operating Temperature	0℃ ~ +40℃
Relative Humidity	30% ~ 90%
Rated Lifting Load	175 kg

Packing & Overall Dimensions

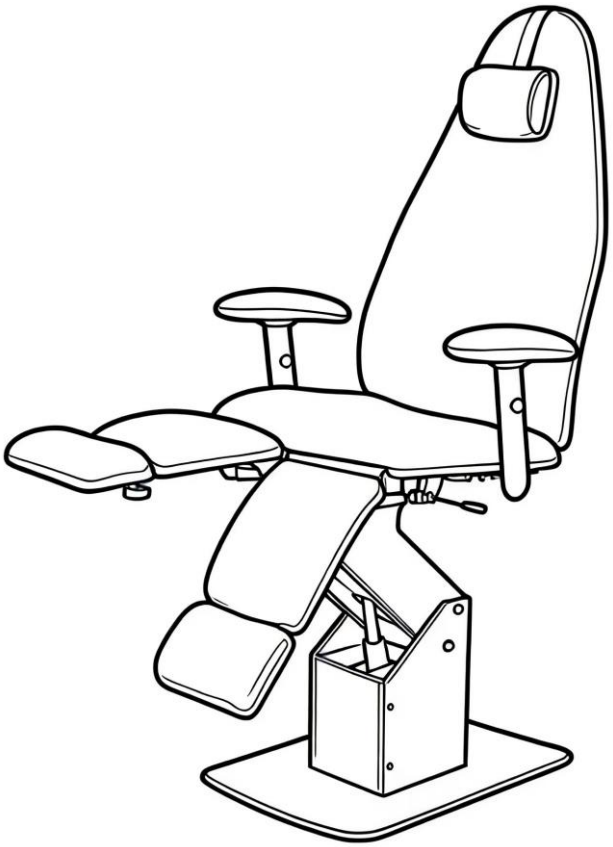
- Packing volume: 0.54 m^3
- Net Weight(NW): 65 kg
- Gross Weight(GW): 76 kg
- Outer carton dimension: 1250 mm × 610 mm × 700 mm
- Bed overall length: 1720 mm; overall width: 765 mm
- Max. backrest tilt angle: 50°
- Max. footrest tilt angle: 70°
- Bed height adjustment range: 630 mm ~ 820 mm

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